


## 防カビ剤の分析

### Thiabendazole, *o*-Phenylphenol, Diphenyl

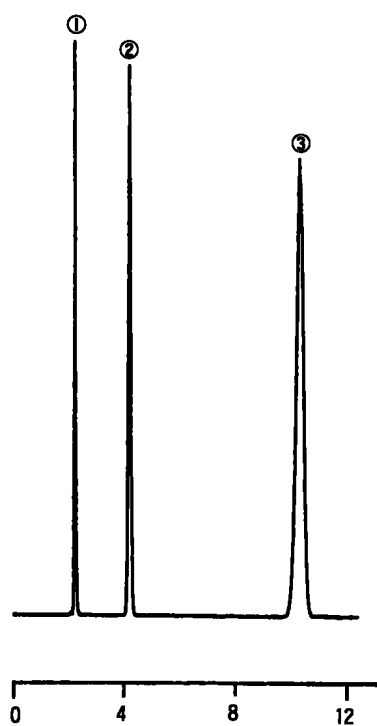

#### Analytical Conditions

Column	: L-column ODS 4.6mm × 150mm
Mobile phase	: 20mM NaH <sub>2</sub> PO <sub>4</sub> /Acetonitrile =40/60
Flow rate	: 1mL/min
Temp.	: 40
Detector	: UV254nm 0.16AUFS
Sample	: Thiabendazole (TBZ) <i>o</i> -Phenylphenol (OPP) Diphenyl (DP)