

## 接着剤中有機溶剤のヘッドスペースガスの分析

## Analysis of Organic solvents in adhesive

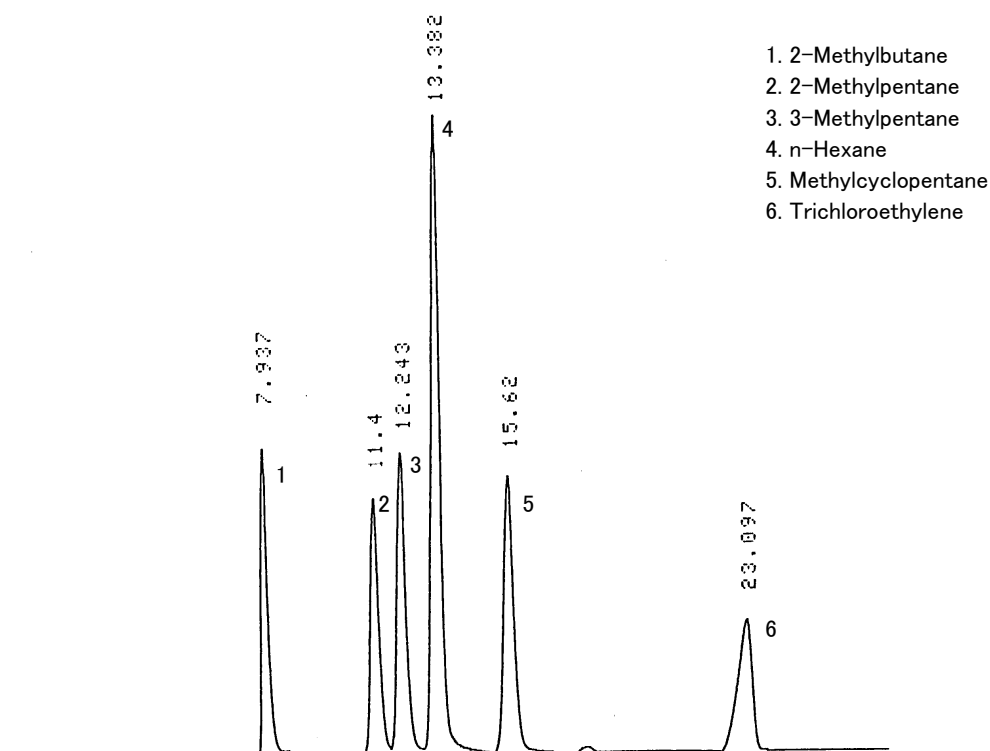


Fig. Chromatogram of Organic solvents in adhesive.

## 【Analytical conditions】

Column: **G-column G-100** 40 m L. × 1.2 mm I.D. 3 μm film thickness  
Flow rate: N<sub>2</sub> 0.2 kg/cm<sup>2</sup>  
C.Temperature: 50°C  
Detection: FID  
Injection volume: 2 mL (Head space gas)  
System: GC-8A (Shimadzu Co.)