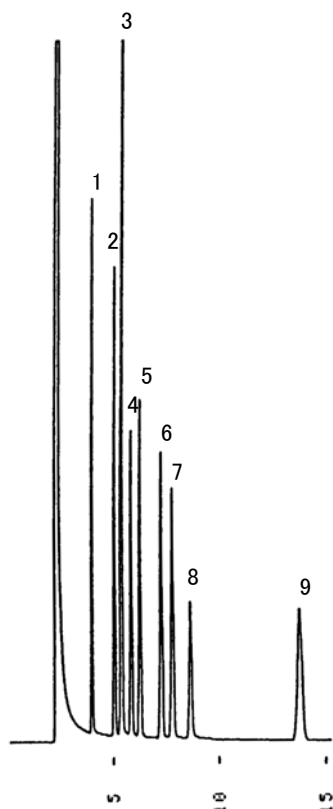


**フェノール類の分析****Analysis of Phenols**

1. Phenol
2. o-Cresol
3. m,p-Cresol
4. Guaiacol (o-Methoxyphenol)
5. 2,6-Dimethylphenol
6. 2,5-Dimethylphenol
7. 3,5-Dimethylphenol
8. p-Chlorophenol
9. p-t-Butylphenol

(100 ppm each)

Fig. Chromatogram of Phenols.

**【Analytical conditions】**

Column: **G-column G-100** 40 m L. × 1.2 mm I.D. 2 μm film thickness  
Flow rate: He 22 mL/min  
C.Temperature: 120°C  
Detection: FID  
Injection volume: 1 μL  
System: -