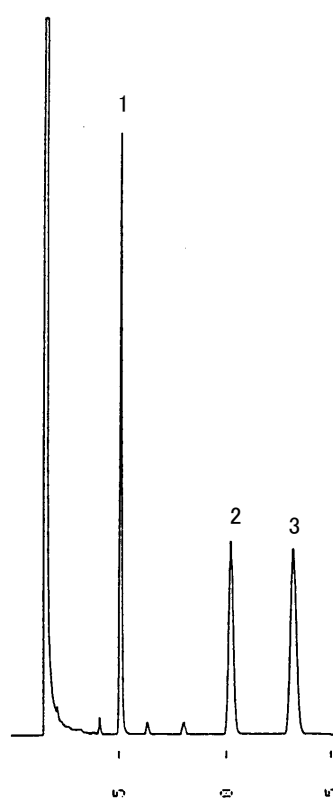


可塑剤の分析**Analysis of Plasticizer**

1. Di(2-ethylhexyl) adipate
2. Di(2-ethylhexyl) azelate
3. Di(2-ethylhexyl) sebacate

(200 ppm each)

Fig. Chromatogram of Plasticizer.

【Analytical conditions】

Column: **G-column G-100** 40 m L. × 1.2 mm I.D. 0.5 μ m film thickness
Flow rate: He 28 mL/min
C.Temperature: 240°C
Detection: FID
Injection volume: 1 μ L
System: -