



アシュラム、オキシン銅の分析

Analysis of Asulam, Oxine copper

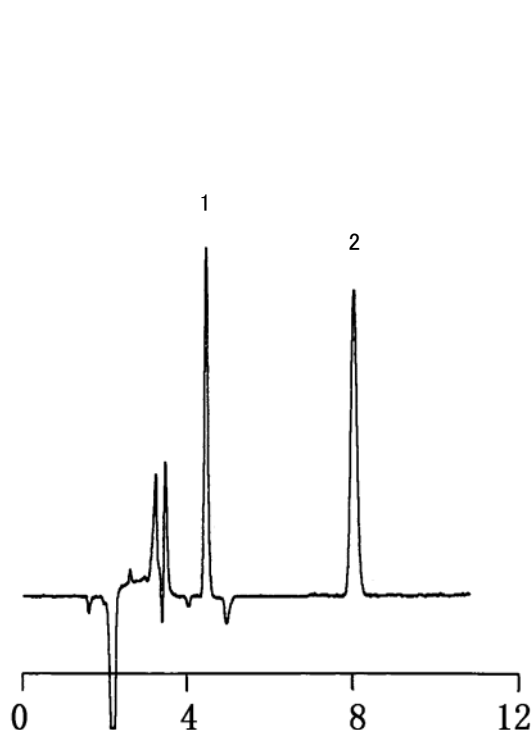

1. Asulam (0.2 mg/L)

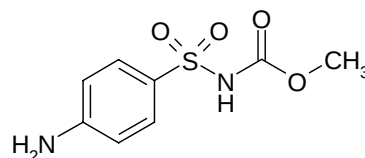

2. Oxine copper (0.04 mg/L)

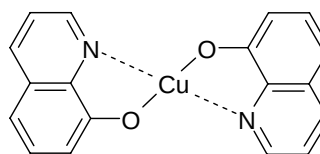


Fig. Chromatogram of Asulam, Oxine copper.

【Analytical conditions】

Column:	L-column ODS (C18, 5 μ m, 12 nm) 4.6 mm I.D. $\times$ 150 mm L.
Mobile phase:	CH ₃ CN/5 mM Sodium octansulfonate + 20 mM H ₃ PO ₄ (18/82)
Flow rate:	1 mL/min
Temperature:	40°C
Detection:	UV 250 nm
Injection volume:	20 μ L
System:	—