


防かび剤
Fungicides

Key words : 防かび剤 ポストハーベスト農薬 食品添加物
Column : USP category: L1


[Analytical conditions]

Column : L-column ODS (C18, 5 μ m, 12 nm), 4.6 mm I.D. \times 150 mm L.; Cat. No. 622070
Eluent : CH₃OH/20 mM NaH₂PO₄ in H₂O (70/30)
Flow rate : 1 mL/min
Temperature : 40°C
Detection : UV 230 nm
Injection volume : 1 μ L
System : LC-9A series (Shimadzu Co.)

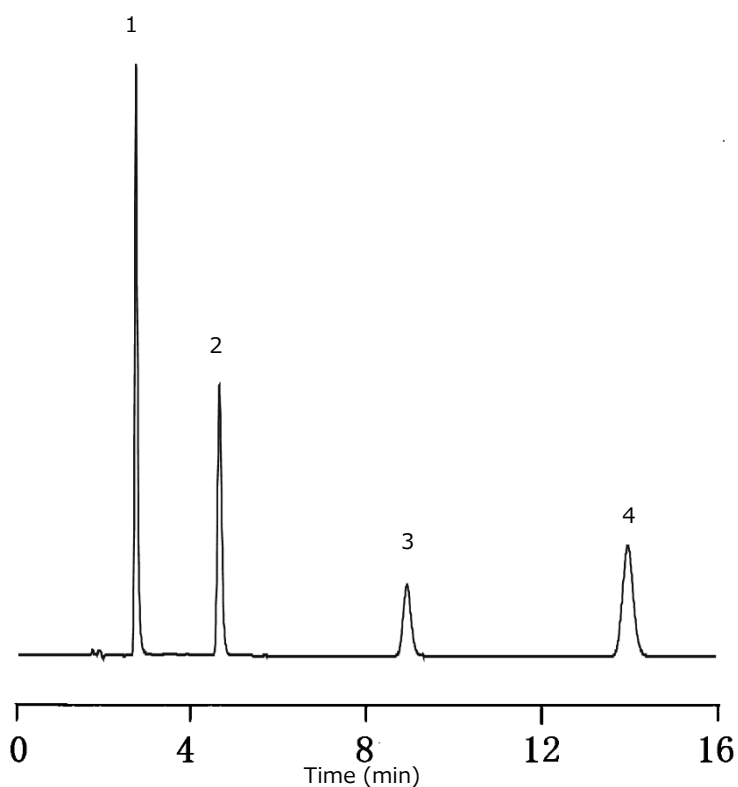


Sample:


1. Thiabendazole (TBZ) (50 mg/L)
チアベンダゾール

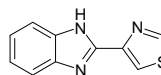


2. o-Phenylphenol (OPP) (50 mg/L)
オルトフェニルフェノール

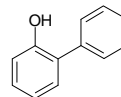

3. Imazalil (50 mg/L)
イマザリル

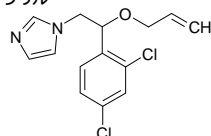

4. Diphenyl (DP) (50 mg/L)
ジフェニル

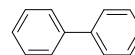

Sample solvent: CH₃CN