


甘味料
Sweeteners

Key words : 甘味料 食品添加物
Column : USP category: L1


[Analytical conditions]

Column : L-column ODS (C18, 5 μ m, 12 nm), 4.6 mm I.D. \times 150 mm L.; Cat. No. 622070
Eluent : CH₃CN/20 mM Phosphate buffer pH 2 (15/85)
Flow rate : 1 mL/min
Temperature : 40°C
Detection : UV 210 nm
Injection volume : 1 μ L
System : LC-9A series (Shimadzu Co.)

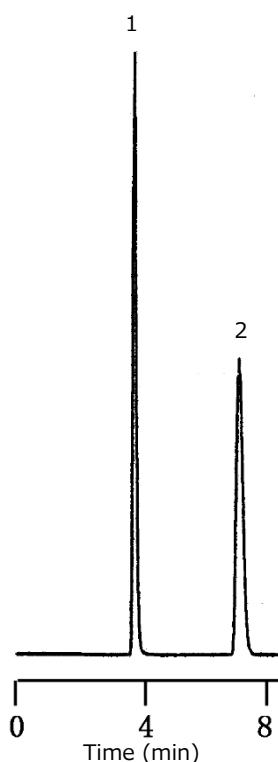

Sample:

1. Saccharin (SA) (200 ppm)
サッカリン

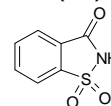

2. Aspartame (APM) (800 ppm)
アスパルテーム

