

Application No. L1032



甘味料 Sweeteners

Key words : 甘味料 食品添加物
Column : USP category: L1

[Analytical conditions]

Column : L-column ODS (C18, 5 μ m, 12 nm), 4.6 mm I.D. $\times$ 150 mm L.; Cat. No. 622070
Eluent : CH₃CN/20 mM Phosphate buffer pH 2.5 (40/60)
Flow rate : 1 mL/min
Temperature : 40°C
Detection : UV 254 nm
Injection volume : 1 μ L
System : LC-9A series (Shimadzu Co.)

Sample:
1. Glycyrrhizic acid (GA) (0.1%)
グリチルリチン酸

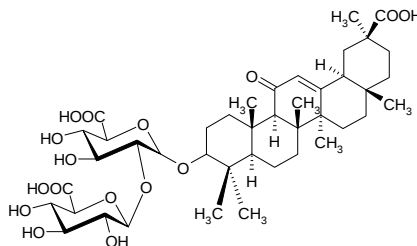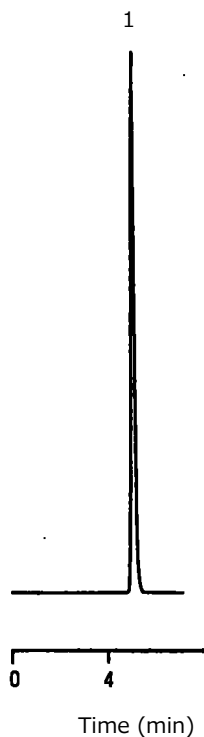

1992 Nishi