


脂肪族アミンの分析

Analysis of Alkyl amines

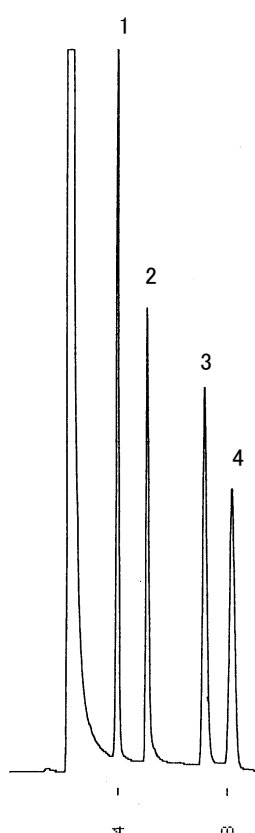

1. Di-n-Butylamine
2. n-Octylamine
3. Di-n-Amylamine
4. Tri-n-Butylamine

(50 ppm each)

Fig. Chromatogram of Alkyl amines.

【Analytical conditions】

Column: **G-column G-205** 40 m L. × 1.2 mm I.D. 5 μm film thickness
Flow rate: He 20 mL/min
C.Temperature: 170°C
Detection: FID
Injection volume: 1 μL
System: -