

## PCB の分析

### Analysis of PCB

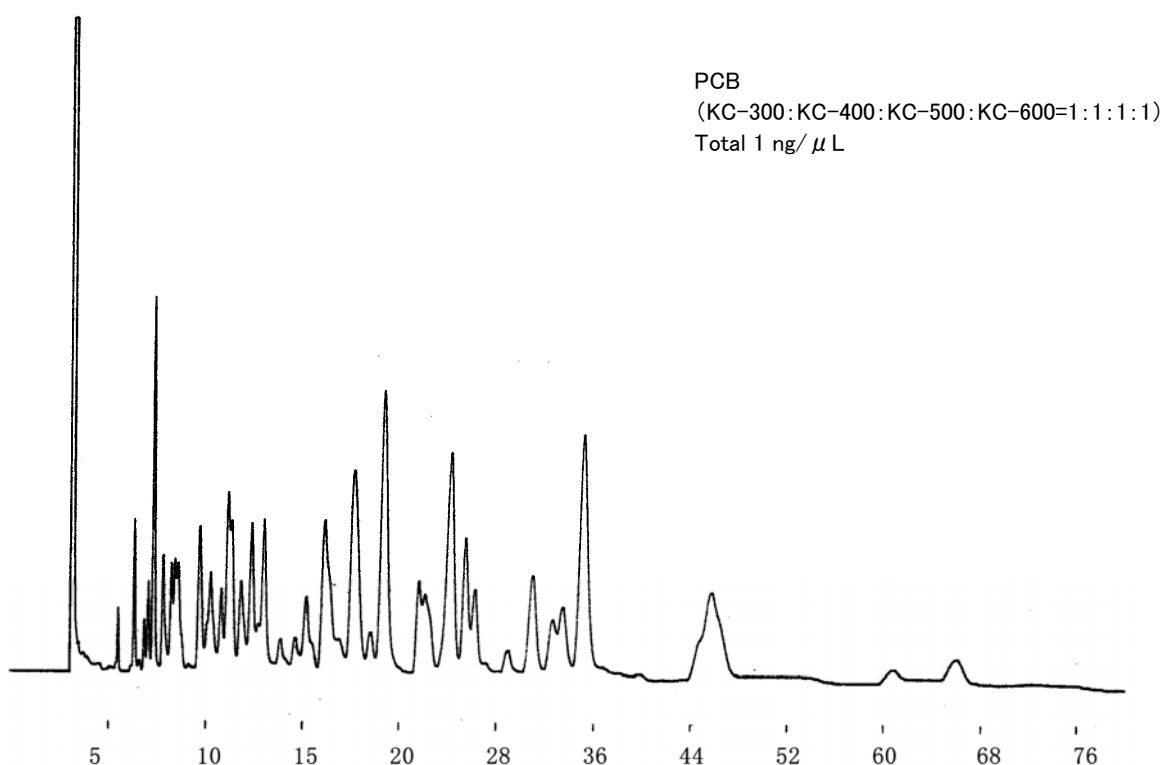


Fig. Chromatogram of PCB.

#### 【Analytical conditions】

Column: **G-column G-250** 40 m L. × 1.2 mm I.D. 0.5 μm film thickness  
Flow rate: He 20 mL/min; Make-up gas: N<sub>2</sub> 30 mL/min  
C.Temperature: 225°C  
Detection: ECD  
Injection volume: 2 μL  
System: -  
\* リードキャピラリー (RN-7, 0.53 mm I.D.) 接続