

有機溶剤の分析

Analysis of Organic solvents

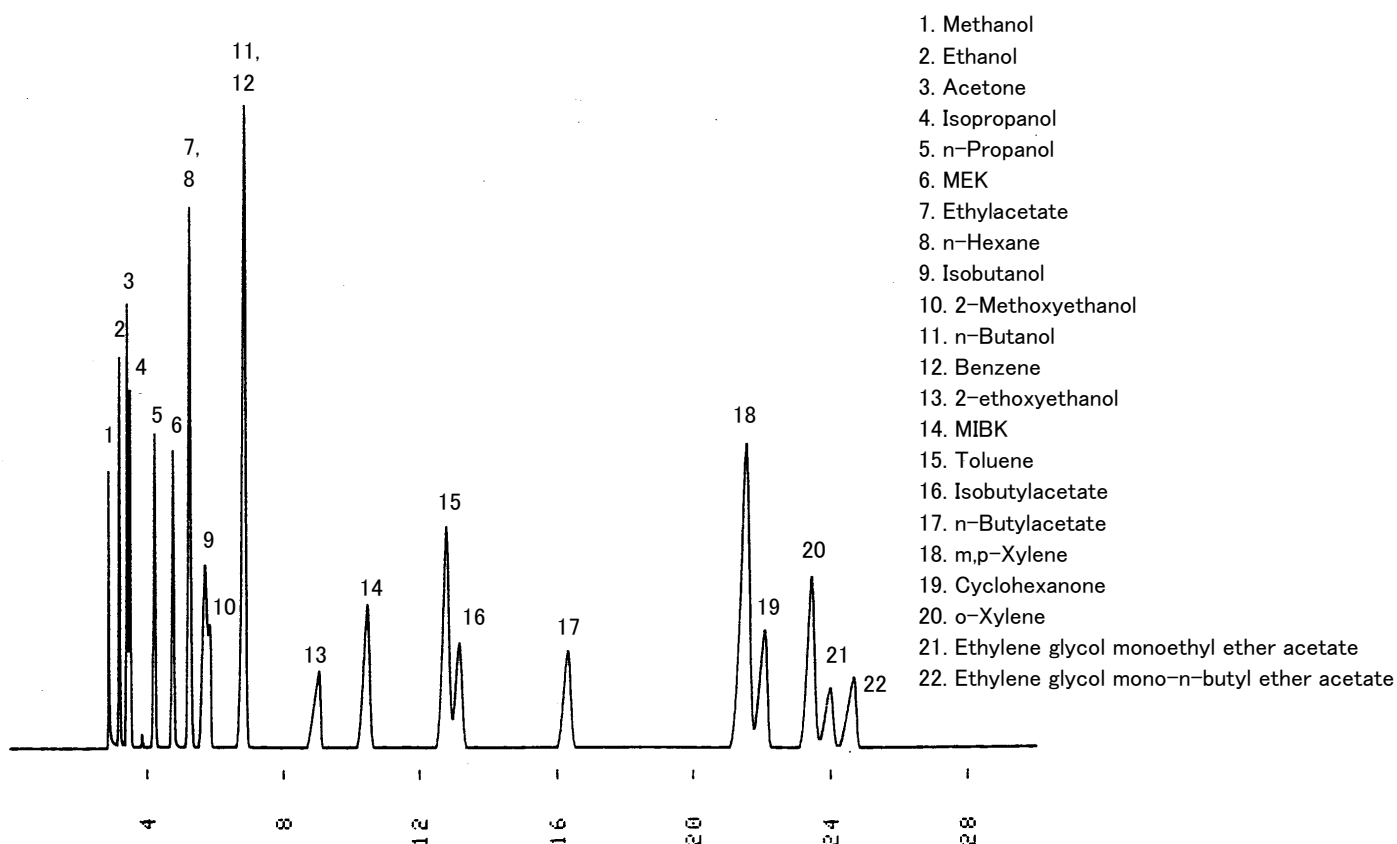


Fig. Chromatogram of Organic solvents.

【Analytical conditions】

| | |
|-------------------|---|
| Column: | G-column G-100 40 m L. × 1.2 mm I.D. 3 μm film thickness |
| Flow rate: | He 17 mL/min |
| C.Temperature: | 50°C (10 min hold) → 2°C/min → 90°C |
| Detection: | FID |
| Injection volume: | 0.5 μL |
| System: | GC-9A (Shimadzu Co.) |